

### ArtOral Modellsystem

More precision – better fit

Enhanced pin system with additional pillar pin.

The big advantage of the ArtOral model system is not only the harmonized combination of materials but also the use of the pillar-pin, which prevents the dies from shifting. As a result a precise position of the die is attained. Additionally, the second pin of the model die system (the soft pin) guarantees an accurate lifting of the plaster stumps.

Item no.	Item	Contents
1406 9999	introductory set	500 piece pillar-pin 1,000 piece smart-pin (incl. 1,000 sleeves) 1 picodent® pinbohrer 1 pin-box blue
1408 0500	pillar-pin, 1.6 mm	500 pieces



### picodent® smart-pin

New, shorter picodent® pin (16.2 mm) with less running resistance. Perfect fit on plaster base.

With even, soft friction. Therefore best suited for the Art-Oral model system.

Item no.	Item	Contents
1406 6001	picodent® smart-pin, 1.6 mm	Set contains 1,000 step pins 1,000 synthetic sleeves
1406 6002	picodent® smart-pin, 1.6 mm	1,000 step pins
1406 6003	picodent® smart-pin, 1.6 mm	1,000 synthetic sleeves




### picodent® pin K "snap"

picodent® step pin-system. A particularly sophisticated pin. With "snap"-effect.

This pin system consisting of cylindrical step pins and synthetic hulls combines the advantages of conical systems with parallel systems.

Item no.	Item	Contents
1406 1000	picodent® pin K "snap", 1.6 mm	Set contains 1,000 step pins 1,000 synthetic sleeves
1406 1002	picodent® pin K "snap", 1.6 mm	1,000 step pins
1406 1001	picodent® pin K "snap", 1.6 mm	1,000 synthetic sleeves

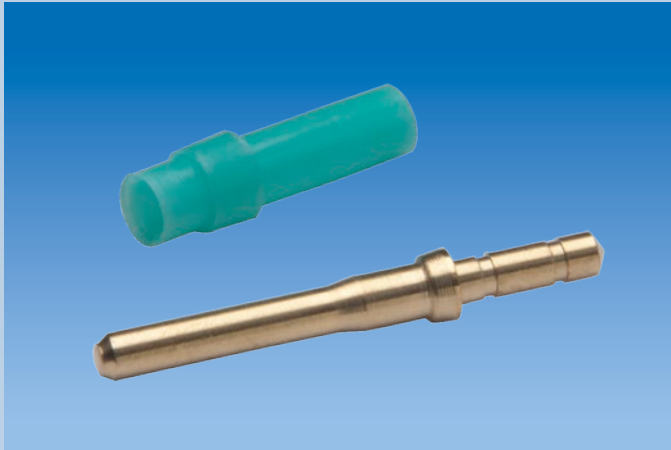

### picodent® soft-pin

Step pins with plastic sleeves and easy fit  
Step pins with soft plastic sleeves, also suitable for Art Oral model system.

Item no.	Item	Contents
1406 2000	picodent® soft-pin, 1.6 mm	Set contains 1,000 step pins 1,000 synthetic sleeves
1406 2002	picodent® soft-pin, 1.6 mm	1,000 step pins
1406 2001	picodent® soft-pin, 1.6 mm	1,000 synthetic sleeves



### picodent® pin MK 2

The “proven” ones with metal sleeve plus guaranteed parallel fit.

The advantage of MK 2 is the collar-like supporting surface of the pins. All pins can thus be fitted in a right angle!

Item no.	Item	Contents
1407 1000	picodent® pin MK 2, 1.6 mm	Set contains 1,000 conical pins 1,000 metal sleeves
1407 1001	picodent® pin MK 2, 1.6 mm	1,000 conical pins
1404 1001	picodent® pin MK 2, 1.6 mm	1,000 metal sleeves



### picodent® pin M

The “proven” ones with metal sleeve  
Combination conical pin + synthetic ring  
or conical pin + metal sleeve.

Item no.	Item	Contents
1404 1000	picodent® pin M, 1.6 mm	Set contains 1,000 conical pins 1,000 metal sleeves
1404 1002	picodent® pin M, 1.6 mm	1,000 conical pins
1404 1001	picodent® pin M, 1.6 mm	1,000 metal sleeves



### picodent® pinbohrer, 1.6 mm

Hard metal drill for ArtOral model system pillar pin, smart pin, pin K, pin MK 2, pin M

For all kinds of drill devices, also available with stronger properties.

Special sizes available on request.

Item no.	Item	Contents
1405 0001	picodent® pinbohrer	1 piece
1405 0002	picodent® pinbohrer mit verstärkter Seele (with stronger properties)	1 piece



### Umkehrformplatte

Opposite matching plate saves time and money

Place the remaining plaster into the opposite matching plate with rubber ring (or magnet with magnet head), place retention and leave to harden.

**The result: a secondary base on stock.**

Item no.	Item	Contents
1600 0016	Umkehrformplatte (8), massive plastic	1 piece




### Formplatten Kunststoff

Coordinated system for producing primary bases

picodent® form plates have been used successfully for years. To produce primary bases, simply add a rubber ring to the plate. Save time and plaster by using form plates that ideally fit the respective size of the model.

Item no.	Item	Contents
1600 0012	Formplatte (7), small, plastic	1 piece
1600 0011	Formplatte (8), standard size, plastic	1 piece

Item no.	Item	Contents
1600 0002	Formplatte (7), small, massive plastic	1 piece
1600 0001	Formplatte (8), standard size, massive plastic	1 piece
1600 0013	Formplatte (9), big, massive plastic	1 piece
1600 0015	Formplatte (11), special size, massive plastic	1 piece

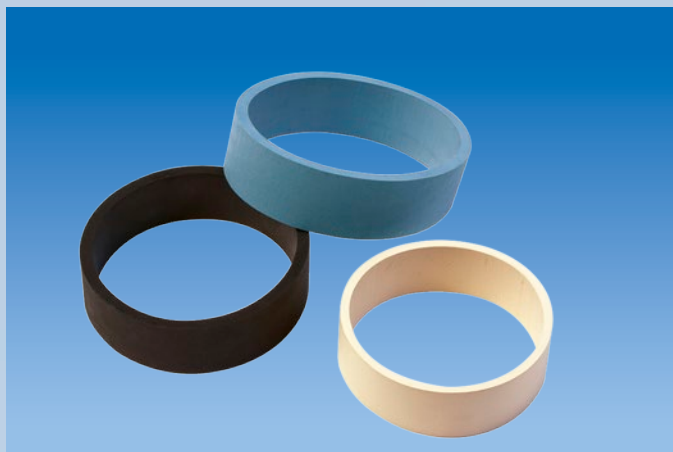

### Gummiringe – for split cast form plates

Matches all pin systems

The standard height (22 mm) of the rubber rings is optimally matched to the plug shaft of our pin systems so that after lifting you can see the top of the pins immediately at the base.

Item no.	Item	Contents
1600 0032	Gummiring (7), small	1 piece
1600 0031	Gummiring (8), standard size	1 piece
1600 0033	Gummiring (9), big	1 piece
1600 0035	Gummiring (11), special size	1 piece

Special height 25 mm for using liquid base stones.




### tip-top sockelformer

For perfect retention of the counterbite model

Produce perfekt antagonist models with the tip-top base former by picodent®. The former consists of high quality rubber base former with integrated retention.

The tip-top base former enables quick work and ensures perfect retention to the antagonist model.

Item no.	Item	Size	Contents
1410 0002	tip-top sockelformer	small	1 piece
1410 0001	tip-top sockelformer	large	1 piece

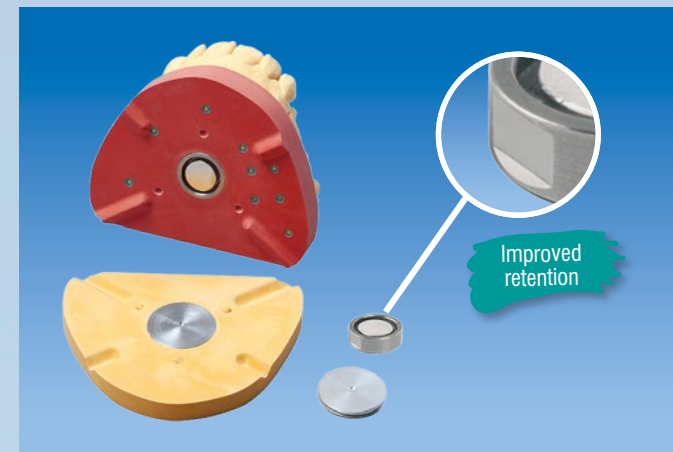

### picolino Magnetsystem

Rust-proof and reusable

Extra flat and small magnet knobs, magnets and catch plates with extreme adhesion for a flat plaster split cast with particularly low tension on the model. Ideal in case of limited space in the articulator. The picolino magnetic system is made of stainless steel and reusable.

Item no.	Item	Contents
1600 3150	picolino Magnete (magnets) Ø 10 mm, thickness 3 mm	50 pieces
1600 3250	picolino Magnettöpfe (magnet knobs) Ø 15 mm, thickness 5 mm	50 pieces
1600 3350	picolino Haftplatten (catch plates) Ø 20 mm, thickness 2 mm	50 pieces
1600 3001	picolino Magnetheber (magnet lifter)	1 piece
1600 3999	picolino Starter-Set	Set consists of: · 20 pcs. picolino Magnete · 20 pcs. picolino Magnettöpfe · 20 pcs. picolino Haftplatten · 1 pc. picolino Magnetheber · 10 pcs. Zentrierkegel






### Magnete · Töpfe · Haftplatten

High adhesive power – reusable. Can also be used with the zeiser® model system.

Item no.	Item	Contents
1600 0042	Magnettöpfe (magnet knobs)	50 pieces
1600 1042	Magnettöpfe Edelstahl, rostfrei (magnet knobs, stainless steel)	50 pieces
1600 0044	Magnete (magnets)	50 pieces
1600 0046	Haftplatten (catch plates)	50 pieces
1600 1046	Haftplatten Edelstahl, rostfrei (catch plates, stainless steel)	50 pieces
1600 0101	Haftplatten Ø 35 mm (catch plates)	100 pieces
1600 0047	Magnetheber (magnet lifter)	1 piece

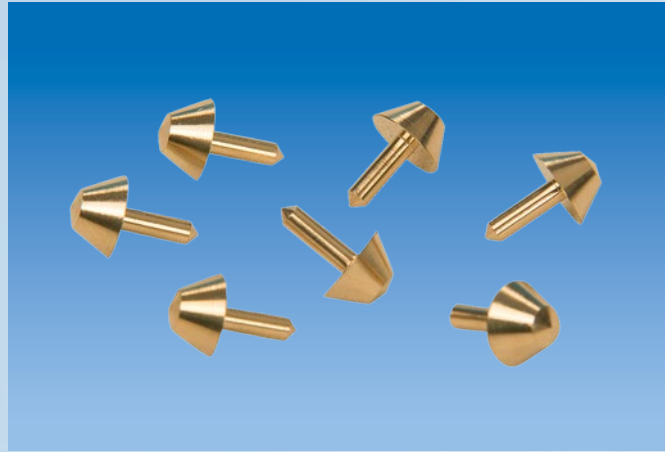

### picodent® zentrierkegel

Provide clean models

Drilling with 1.6 mm pin-drill protects against particles of dirt on the underside of the model.

Item no.	Item	Contents
1600 0051	picodent® zentrierkegel	50 pieces



### picodent® picocryl

Cyanacrylate glue

Low viscosity and very fast. Ideally suited as a pin glue.

Item no.	Item	Contents
1505 0035	picodent® picocryl	10 ml bottle