



WITH NEW FORMULA



WITH NEW FORMULA



pico-pattern

Self-curing universal modelling resin

- Casting of telescopic and conical crowns
- Fabrication of resin stumps
- Casting of individual attachments and bars
- Fixation of splints e.g. of soldered models
- Casting of inlays, onlays, retentions and clasps
- Casting of superstructures

Item no.	Item	Contents
1825 0002	pico-pattern Set	100 g powder 100 ml liquid 1 x pico-pattern brush 2 x measuring cup 5 ml
1825 0001	pico-pattern Metallständer	1 piece
1825 0010	pico-pattern Pulver	100 g bottle
1825 0050	pico-pattern Pulver	500 g bottle
1825 0100	pico-pattern Liquid	100 ml bottle
1825 0250	pico-pattern Liquid	250 ml bottle
1825 0005	Silikonanmischbecher	2 pieces á 5 ml
1825 0006	pico-pattern-Pinsel	10 pieces

picovest® royal NF

Graphite-free phosphate-bound investment material for the crown and bridge technique.

Especially suited for NEM and optimised for rapid casting.

Mixing ratio for 100 g powder:	24 ml liquid
Mixing ratio for 150 g powder:	36 ml liquid
Stirring time under vacuum:	60 seconds
Working time:	approx. 5 minutes
Setting expansion:	approx. 2,9 %

Item no.	Item	Contents
1206 1500	picovest® royal NF Starterbox	12 bags 150 g each 500 ml liquid 50 ml measuring cylinder
1206 3015	picovest® royal NF	30 bags 150 g each
1206 0225	picovest® royal NF	5 kg bag
1206 1010	picovest® royal NF Liquid	1.000 ml bottle

picovest® universal NF


Phosphate-bound precision investment material.

For the crown and bridge technique as well as for all pressable ceramics and optimised for rapid casting

Mixing ratio for 100 g powder:	25 ml liquid
Mixing ratio for 150 g powder:	38 ml liquid
Stirring time under vacuum:	60 seconds
Working time:	approx. 5 minutes
Setting expansion:	approx. 2,0 %

Item no.	Item	Contents
1205 1500	picovest® universal NF Starterbox	12 bags 150 g each 500 ml liquid 50 ml measuring cylinder
1205 0025	picovest® universal NF	5 kg bag
1205 2550	picovest® universal NF	25 bags 150 g each
1205 5015	picovest® universal NF	50 bags 150 g each drum with carrying handle
1205 1010	picovest® universal NF Liquid (liquid)	1.000 ml bottle
1205 1020	Spezial-Liquid für picovest® universal NF (Special Liquid - pink)	1.000 ml bottle

NEW

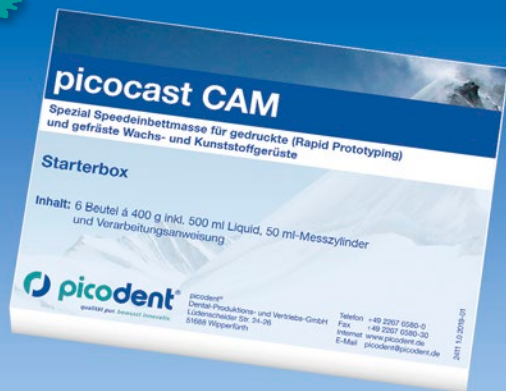

WITH NEW FORMULA



picocast CAM

Special speed investment material for printed (rapid prototyping) and milled wax and resin frameworks

It is optimised for the speed casting techniques and is suitable for casting all types of dental alloys (except Titanium).

Mixing ratio for 100 g powder:	20 ml liquid
Mixing ratio for 400 g powder:	80 ml liquid
Stirring time under vacuum:	60 seconds
Working time:	approx. 5 minutes
Max. curing expansion:	approx. 3 %

Item no.	Item	Contents
1211 6400	picocast CAM Starterbox	6 bags 400 g each 500 ml liquid 50 ml measuring cylinder
1211 5040	picocast CAM	50 bags 400 g each
1211 1000	picocast CAM Liquid	1.000 ml bottle

picocast speed NF

Investment material for model casting technique.

The rapid filling investment material for model casting work. Optimised for speed-casting.

Mixing ratio for 100 g powder:	20 ml liquid
Mixing ratio for 400 g powder:	80 ml liquid
Stirring time under vacuum:	60 seconds
Working time:	approx. 5 minutes
Setting expansion:	approx. 0,8 %

Item no.	Item	Contents
1212 1500	picocast speed NF Starterbox	12 bags 200 g each 500 ml liquid 50 ml measuring cylinder
1212 0450	picocast speed NF	4.5 kg bucket
1212 1845	picocast speed NF	4 x 4.5 kg bucket
1212 5040	picocast speed NF	50 bags 400 g each
1212 1200	picocast speed NF	100 bags 200 g each
1212 1010	picocast speed NF, Liquid	1.000 ml bottle

picodent® easysoft

Gypsum-based investment material

For hot and cold ovens

Special investment material for all precious alloys. Perfect no matter whether the investment ring with fleece inlay is placed in a cold or hot oven.

Mixing ratio for 100 g powder:	32-37 ml distilled water
Mixing ratio for 130 g powder:	42-48 ml distilled water
Stirring time under vacuum:	60 seconds
Working time:	approx. 7 minutes
Final set after:	approx. 60 minutes
Setting expansion:	0.4-0.8 % depending on liquid concentration
Thermal expansion:	1.2 %

Item no.	Item	Contents
1200 0006	picodent® easysoft	6 kg bag




picodent® easyvlies

Asbestos-free investment ring

Strong interior lining for casting muffles.

Item no.	Item	Contents
1203 0055	picodent® easyvlies	roll of 55 mm x 1 mm x 25 m

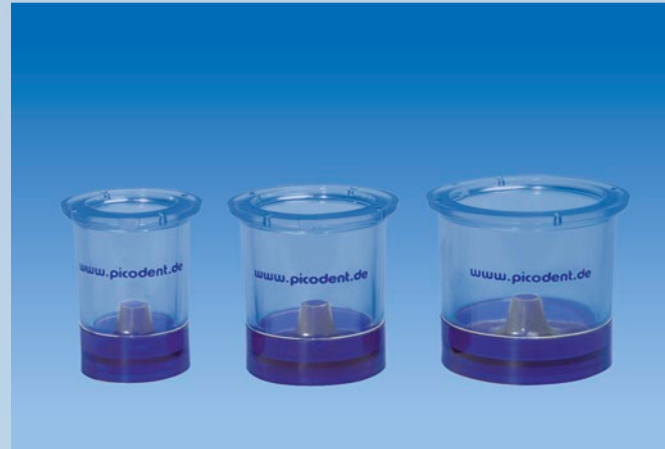

picodent® „Speed-Muffle“

Saving with system

Metal ring-free work. Three dimensional expansion. Perfect fit.

Item no.	Item
1203 0002	picodent® „Speed-Muffle“ - Größe 2
1203 0003	picodent® „Speed-Muffle“ - Größe 3
1203 0006	picodent® „Speed-Muffle“ - Größe 6